

# Qufora® Leg bag, non-sterile


## Discreet. Comfortable. Easy-to-use.

The **Qufora®** Leg bags are designed to provide discreet storage of urine under clothing and can be attached to either a urinary catheter or an urisheath. The **Qufora®** Leg bags are mainly used during the day and are easily fixed to the leg using straps.

### Product benefits

- The **Qufora®** Leg bag has a universal conical inlet connector making it easy to attach a urinary catheter or urisheath to the bag.
- A non-return valve makes sure that urine will not flow back into the tube and infection risks are thereby highly reduced.
- The multi-chamber system ensures that the **Qufora®** Leg leg bag stays discreet even when it's full.
- Soft backing, modern design and shape add a high degree of comfort for the end-user.
- The patented Q-Flow valve is easy and safe to operate also for people with reduced finger dexterity. Turning it 90° the valve is fully open. At 180° the valve is still fully open and now placed in a safe position without any risk of closing unexpectedly during the night.
- The Q-Flow valve offers a silicone tube for easy connection to a urine bag for the night, decreasing risk of infection due to change of urine bag.
- The **Qufora®** Leg bag is available with an adjustable tube, making it possible to personalise the tube's length to fit perfectly for all needs.
- The **Qufora®** Velcro strap system is designed for the **Qufora®** Leg bag System. The straps are adjustable for any individual's needs and the **Qufora®** Leg bag includes eyelets for the velcro straps.

Qufora® Leg bag

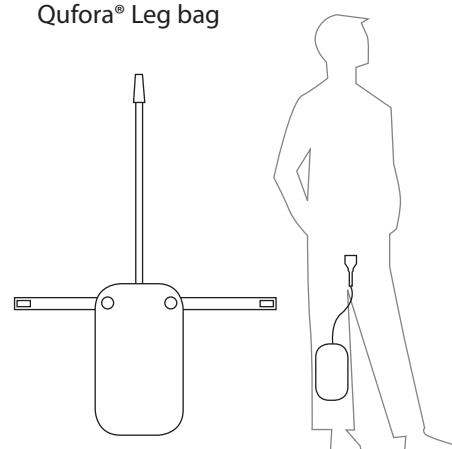

Connector



Q-Flow valve



Scale, multi-chamber, soft backing



# Qufora® Leg bag, non-sterile

## Technical data

### Product description

**Qufora®** Leg bags are available in several versions depending on features like volume size or length of the inlet tube.

The **Qufora®** Leg bags are single-use medical devices.

**All materials are 100% latex- and DEHP-free.**

<b>Qufora® Leg bag features</b>	8 mm diameter inlet tube with a frosted surface Connector cap Transparent foil/white 0,2 mm Non-return valve Scale on the front Non-woven soft backing Eyelets for straps Outlet valve with soft silicone tube	
<b>Material</b>	Connector cap Connector Inlet & Outlet tube Bag body Non-return valve Backing Outlet valve Valve tube	PE PVC (DEHP-free) PVC (DEHP-free) PVC (DEHP-free) PVC (DEHP-free) PET PVC PP + PE
<b>Sterilisation</b>	Non sterile	
<b>Logistics</b>	Product	Single packed
Packaging	Retail box	H: 120mm X W: 160mm X D: 230mm
	Shipper	H: 240mm X W: 330mm X D: 615 mm
Shelf life	IFU	Included in each retailbox
		5 years
Storage		Dry in temperatures between 0°C and +35°C, protected from direct sunlight.
		Not to be exposed to heavy pressure
<b>GMDN</b>	Non sterile	58924
<b>Product classification</b>	Non sterile	Medical devices acc. to MDD 93/42/EEC, CE class I
<b>Applied standard</b>	ISO 13485 ISO 8669-2	Quality Management System Urine collection bags — Part 2: Requirements and test methods

