


Quality time for better care

"Quality time for better care" is the brand promise of Terumo Medical Care Solutions

HOSPITAL CARE

# Surflo™

## Winged Infusion Set with Protector Product Information

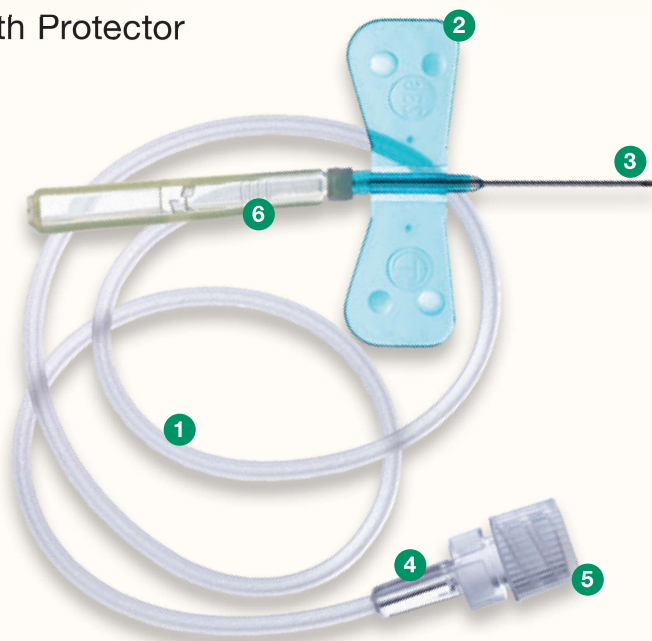

The Surshield protector is designed to allow the user's fingers to remain behind the needle point and thus minimizing the risk of an accidental needle-stick injury. Once the protector shield is activated, the needle is permanently covered, allowing the end-user to safely dispose of it.

Description	Raw Materials
1 Tube	Polyvinyl chloride
2 Wing hub	Polyvinyl chloride
3 Cannula	Stainless steel
4 Adapter	PMMA
5 Adapter cap	Polypropylene
6 Surshield protector	Polypropylene
Safety hinge	ABS
Protector	Polyethylene

Does not contain components made of DEHP or natural rubber latex.

### Sterilization

With ethylene oxide to a Sterility Assurance Level of at least 10<sup>-6</sup>

### Storage

- Shelf life of 3 years
- The batch number and expiry date are printed on the unit pack
- Do not store at extreme temperature and humidity

Easy activation  
of the safety feature

## The comfort of your patients in your hands

### Product Range

Product Code	Colour	External Cannula Ø		Cannula Length	Dead Space Volume	Tubing Length
SV*S18NL30	Pink	18G	1.20 mm	19 mm	0.45 ml	300 mm
SV*S19NL30	Cream	19G	1.10 mm	19 mm	0.45 ml	300 mm
SV*S21NL09	Green	21G	0.80 mm	19 mm	0.20 ml	90 mm
SV*S21NL30	Green	21G	0.80 mm	19 mm	0.41 ml	300 mm
SV*S22NL09	Black	22G	0.70 mm	19 mm	0.20 ml	90 mm
SV*S22NL30	Black	22G	0.70 mm	19 mm	0.41 ml	300 mm
SV*S23NL09	Blue	23G	0.60 mm	19 mm	0.20 ml	90 mm
SV*S23NL30	Blue	23G	0.60 mm	19 mm	0.40 ml	300 mm
SV*S25NL30	Orange	25G	0.50 mm	19 mm	0.40 ml	300 mm



### Regulatory Status

CE Certification – Class IIa medical device

**Notified Body**  
TÜV SÜD 0123

**Manufacturer**  
TERUMO MEDICAL PRODUCTS HANGZHOU, M4-9-5,  
Economic & Technological Development Zone, Hangzhou, 310018 China

**EC Representative**  
TERUMO EUROPE N.V., Interleuvenlaan 40, 3001 Leuven, Belgium

### Precautions

- Do not use if unit pack is damaged
- For single use only
- Dispose of safely after single use to avoid risk of infection

### Packaging

Unit Box	Shipping Carton		
	Quantity	Gross Weight	Dimensions (L x W x H)
50	500	3,2 kg	510 x 257 x 260 mm
50	500	3,2 kg	510 x 257 x 260 mm
50	500	2,9 kg	510 x 257 x 260 mm
50	500	3,1 kg	510 x 257 x 260 mm
50	500	2,9 kg	510 x 257 x 260 mm
50	500	3,1 kg	510 x 257 x 260 mm
50	500	2,9 kg	510 x 257 x 260 mm
50	500	3,1 kg	510 x 257 x 260 mm
50	500	3,1 kg	510 x 257 x 260 mm

Single unit: Blister lid paper / thermoformable blister film

### Choose from a wide range from 18G up to 30G

The outside diameter of the needle can be recognized by the ISO colour-coded hub.

Your Terumo partner



**TERUMO EUROPE N.V.**

Terumo Corporation    Terumo Europe N.V.  
+81 3 6742 8500    +32 16 38 12 11

[www.terumo-europe.com](http://www.terumo-europe.com)

All brand names are trademarks or registered trademarks of TERUMO CORPORATION and their respective owners. Published by Terumo Europe N.V.