




Floswitch®
Sterile, single use


- Floswitch® acts as an ON/OFF device on infusion lines, stopping the flow in order to carry out line procedures in closed circuit.
- Designed to reduce risk of patient emboli and contamination of hospital staff to blood.
- The Floswitch® is designed to increase safety.
- Floswitch® can be connected to any catheter fitted with a Luer connector.
- Winged to facilitate attachment.

Floswitch® closed

External compression of the extension tubing by means of a steel ball, provides risk-free unambiguous opening and closing of the infusion line.

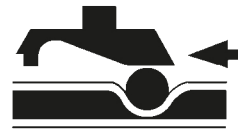

Floswitch® open

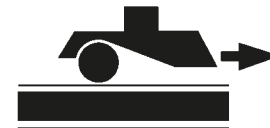

ORDERING INFORMATION

PRODUCT CODE	DESCRIPTION	BOX/CASE CONTENT (UNITS)
681033	Floswitch®	25/250

Latex free

CE marked

Class IIa according to the European Medical Devices Directive 93/42/EEC as amended by Directive 2007/47/EC Shelf life: 5 years

Sterilisation method: Ethylene Oxide