#### COMON



Intuitive operation interface: Configured with 3"LCD



High adaptability design: Multifunctional connector for RS232, AC, nurse call.

Sleek ergonomic design:

• Longer Battery Life: No less than 12hrs\*

\*Under specific conditions.









#### Shenzhen Comen Medical Instruments Co., Ltd.

Add: Floor 10, Floor 11 and Section C of Floor 12 of Building 1A &Floor 1 to Floor 5 of Building 2, FIYTA Timepiece Building, Nanhuan Avenue, Matian Sub-district, Guangming District, Shenzhen, Guangdong, 518106,P.R. China Tel:+86-755-2640 8879 Fax:+86-755-2643 1232 Website: en.comen.com E-mail: info@szcomen.com

P/N: EN-M260-V-8P-20220621-V1.0



M260-V Infusion Pump

### Our Latest Frugal Innovations

Our next generation of economical pumps goes beyond cost.



**EN1789** certified for use during transport and E&R scenarios. Shielded from harsh environments with a validated **IP44** rating.





# What makes it new? It's the way you use it.



#### **Upgraded:**


Bigger screen. Lighter Weight

>>>> Wider flow rate ranges

Higher sensitivity and accuracy

Larger storage capacity

Longer battery endurance

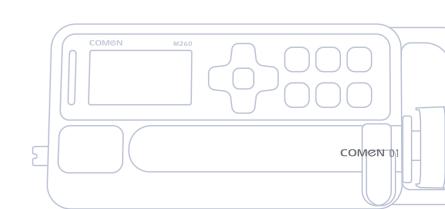


## The New Multifunction Comprehensive Solution


#### Functionally Extensive. Clinically Inclusive.


10 different working modes with wider setting ranges:

Rate/Time/Weight/Drip<sup>[1]</sup>/Dose Time/Intermittent/Ramp/Micro/First Dose/Sequential

- Safer for pediatric/neonatal use
- Smooth & flexible workflow up to 9 phases of ramp up/down and up to 10 sequential setups available
- Rates of up to **2000**ml/h
- Infusion Time up to  $\sim 100$  hrs<sup>[2]</sup>

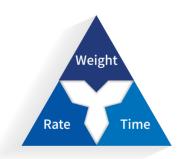


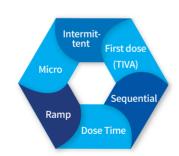

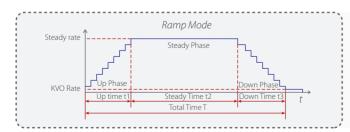

Recommended use for some vasodilators, (eg. Nitroglycerin and triclopidine), to prevent abrupt discontinuation causes rebound vasoconstriction and angina attacks strains.

[2] Time set up 00:00:01~99:59:59 hh:mm:ss

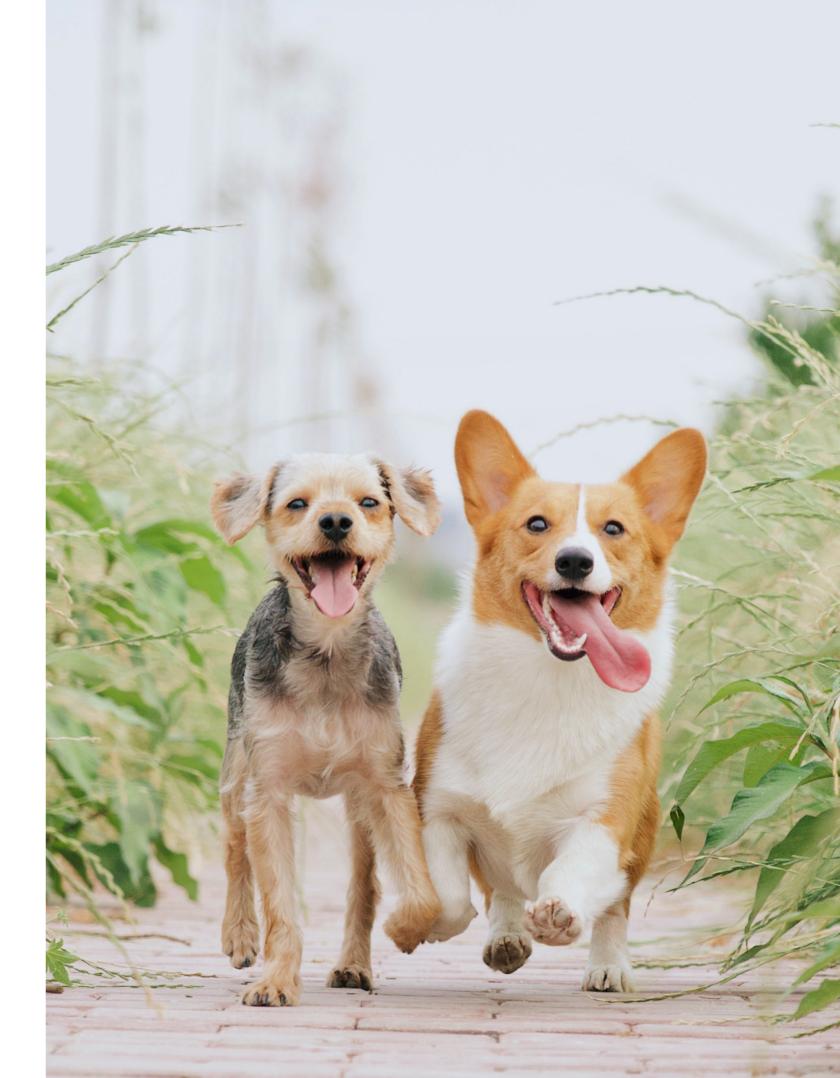

<sup>\*[1]</sup> Drip mode only applicable for ME660-V.

### M260-V

#### **Enhanced Sensitivity and Accuracy** with Strong Functionality.


0.01 Accurate increments

±1.8% Infusion accuracy<sup>[3]</sup>




0.5ml Maximum bolus under

#### Robust emergency support with excellent portability.

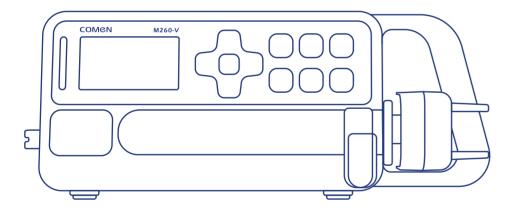

- Super lightweight: only **1.6**kg<sup>[4]</sup> with handle
- Longer battery endurance: up to 12hrs[5]
- Adjustment for multi-environment: 10 toggleable levels for screen brightness and alarm volume
- QuickStart protocol to boost aid-giving
- \*[3] When infusion rate≥1.00ml/h
- [4] Including built-in battery
- [5] With optional 5100mAh battery, under rates of 5ml/h

#### SafeGuard: **Risk Mitigation Solutions.**

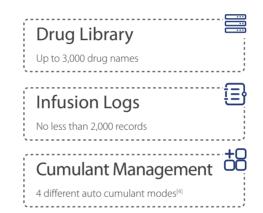

#### **Patented DPS system**

Up to 15 selectable pressure release triggers

Dynamic Pressure Monitoring and Automatic Pressure Release system.

Real-time dynamic pressure detection





\*[6] Auto cumulant modes include: Recent cumulant/Timed cumulant/Interval cumulant/24h cumulant.