



# EN-V5 Vet Infusion Pump



## with Heater

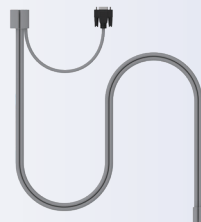


It can only be used with the  
**EN-2203 Heater of ENMIND**

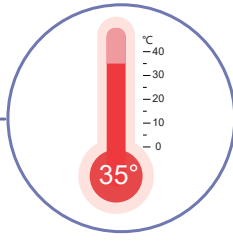



**Flexible Adjustment:**  
30~40°C continuously adjustable, temperature control accuracy  $\pm 1.0^{\circ}\text{C}$ .



**Multiple Protections:**  
The heating will be automatically turned off when the temperature exceeds 50°C to ensure the safety of infusion.

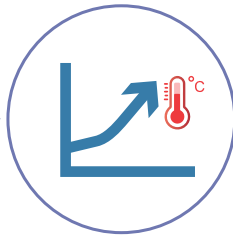


**Rapid Heating:**  
4 sets of heating units quickly reach the set temperature.

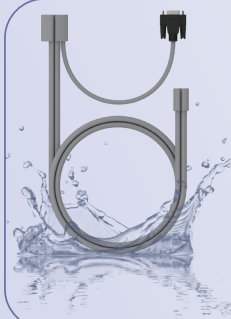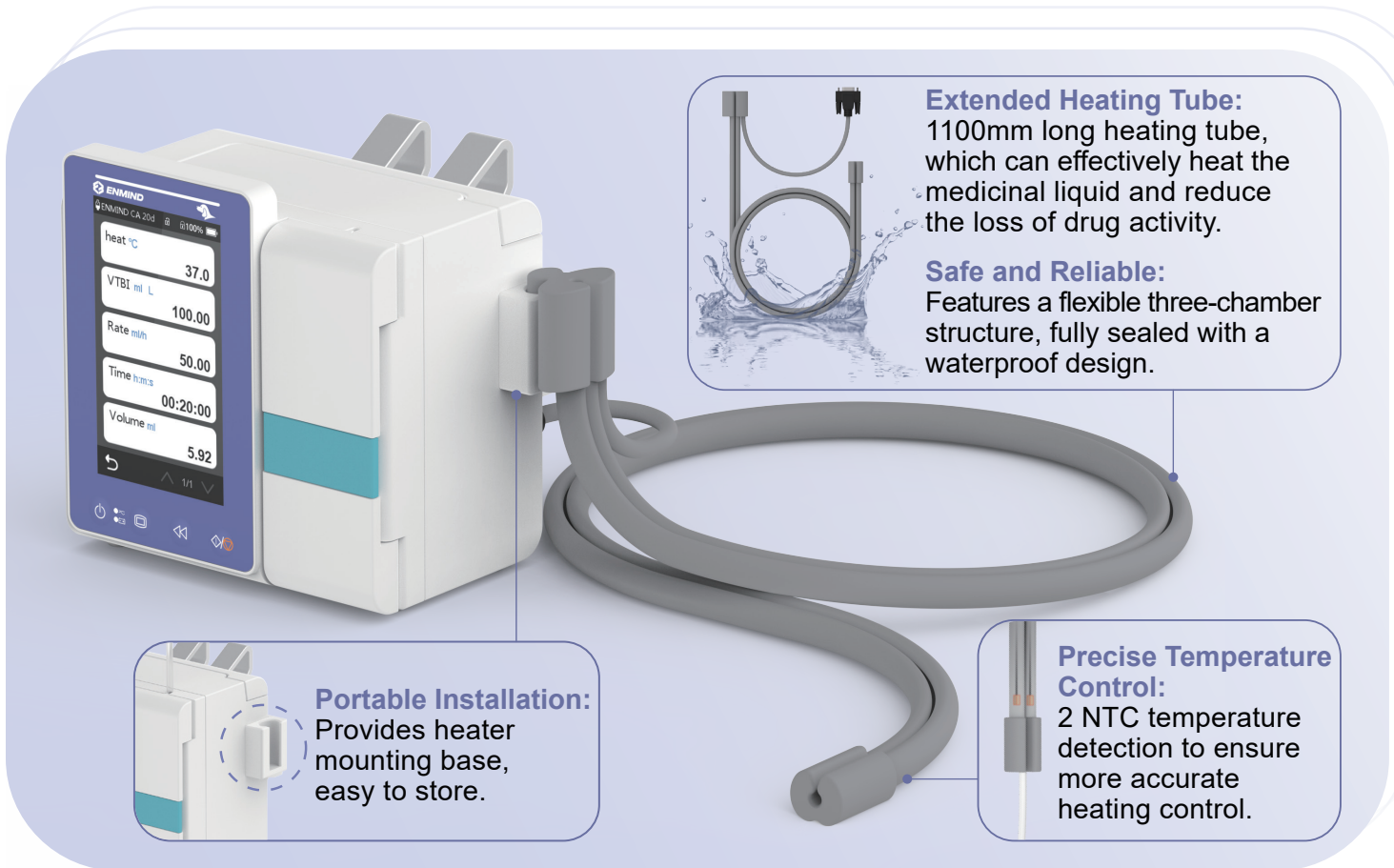

**Extended Heating Tube:**  
1100mm long heating tube, which can effectively heat the medicinal liquid and reduce the loss of drug activity.

**Safe and Reliable:**  
Features a flexible three-chamber structure, fully sealed with a waterproof design.

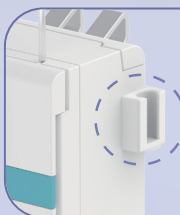

**Portable Installation:**  
Provides heater mounting base, easy to store.



**Precise Temperature Control:**  
2 NTC temperature detection to ensure more accurate heating control.

