


HOSPITAL CARE

Versatus - WP

Winged and Ported I.V. Catheter
Product Information

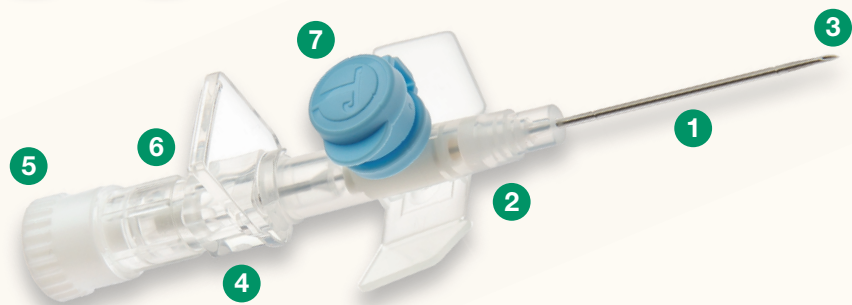

- The flashback chamber visually confirms that the catheter is correctly inserted in the vein
- The sharp needle tip minimizes trauma to the vessel during use
- Wide range from 14G up to 26G

Description	Raw Materials
1 Catheter tube	Polyurethane
2 Winged catheter hub Caulking pin	Polypropylene Polyoxymethylene
3 Cannula (internal needle)	Stainless steel
4 Needle hub	Styrene butadiene copolymer
5 Luer-lock cap	Polypropylene
6 Vent fitting filter Vent fitting body	Polyethylene Styrene butadiene copolymer
7 Port Closure Injection port valve	Polypropylene + master batch Silicone rubber
Protector	Polyethylene
Unit pack	Adhesive coated lid paper. Thermoformed blister film: laminated polyethylene/ polypropylene complex

Does not contain components made of PVC or natural rubber latex.

Designed and
manufactured by skilled
technicians in Europe

Your catheters of choice for a wide range of I.V. applications

Product Range

Product Code	Colour	External Catheter Ø		Internal Catheter Ø	Catheter Length	Flow Rate ml/min.
SR+DU2619PD	Violet	26G	0.64 mm	0.45 mm	19 mm	19
SR+DU2419PD	Yellow	24G	0.74 mm	0.55 mm	19 mm	29
SR+DU2225PD	Blue	22G	0.90 mm	0.65 mm	25 mm	42
SR+DU2032PD	Pink	20G	1.00 mm	0.75 mm	32 mm	59
SR+DU1832PD	Green	18G	1.30 mm	0.95 mm	32 mm	103
SR+DU1845PD	Green	18G	1.30 mm	0.95 mm	45 mm	96
SR+DU1745PD	White	17G	1.50 mm	1.15 mm	45 mm	155
SR+DU1645PD	Grey	16G	1.75 mm	1.35 mm	45 mm	225
SR+DU1445PD	Orange	14G	2.00 mm	1.55 mm	45 mm	290



Packaging

Single Unit Blister pack

Unit Box 50 pieces

Shipping Carton 500 pieces, Gross Weight ± 4.5 kg
Dimensions ± 385 x 379 x 170 mm

Regulatory Status

CE Certification – Class IIa medical device

Notified body

TÜV SÜD Product Service GmbH (0123)

Manufacturer

Delta Med S.p.A. (Italy)

Distributor

TERUMO EUROPE N.V., Interleuvenlaan 40, 3001 Leuven, Belgium

Sterilization

With ethylene oxide to a Sterility Assurance Level of at least 10⁻⁶

Storage

- Shelf life of 5 years
- The batch number and expiry date are printed on the paper of the unit pack
- Keep away from direct sunlight and heat.
- Keep dry
- Fragile, handle with care

Precautions

- Do not use if the package is damaged and consult instructions for use
- Do not re-use / do not re-sterilize
- Do not attempt to re-insert a partially or completely withdrawn needle
- Dispose of safely after single use to avoid risk of infection
- Application time should be consistent with accepted standards of practice

Always consult the product label and the instructions for use for a complete overview of warnings, cautions and/or precautions prior to actual use.

This range complements the wider Terumo product offering including a comprehensive portfolio of both safety and conventional I.V. products which can be found at our web page on www.terumo-europe.com



Your Terumo partner

TERUMO EUROPE N.V.

Terumo Corporation Terumo Europe N.V.
+81 3 6742 8500 +32 16 38 12 11

www.terumo-europe.com

All brand names are trademarks or registered trademarks owned by TERUMO CORPORATION, its affiliates, or unrelated third parties. Published by Terumo Europe N.V.

Image for visualisation purpose only