




## Now with Point-of-Use Dispensing Options

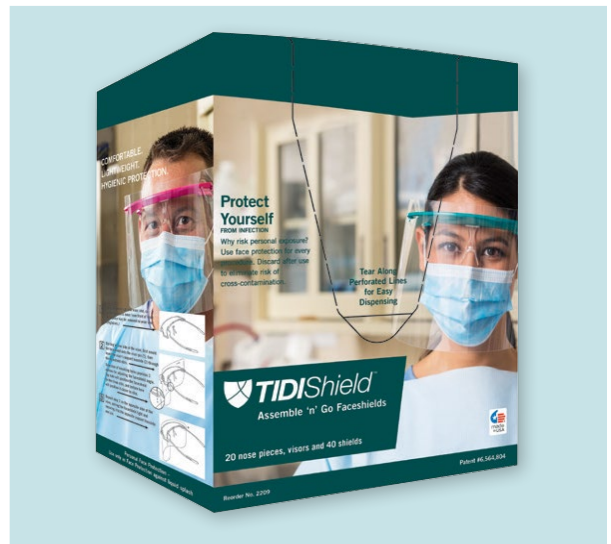

*TIDIShield® Assemble 'n Go™ Face Shields  
PROTECTION and COMFORT*

## PROTECTION IS ESSENTIAL – TIDISHIELD® ASSEMBLE 'N GO™ FACE SHIELDS ARE THE ANSWER

Healthcare professionals work under hazardous and ever-changing conditions. What healthcare workers need is superior face protection from the splash and splatter of blood-borne pathogens and other potentially infectious materials. TIDIShield® Assemble 'n Go™ Face Shields are the answer.

TIDIShield® Assemble 'n Go™ Face Shields are a cost-effective face protection product. Face shields are lightweight, cool and offer comfortable protection that can be worn over prescription eyewear and facemasks.

### THE NEED FOR PROTECTIVE FACEWEAR

- The mucous membranes of the mouth, nose, and eyes are susceptible portals of entry for infectious agents, as can be other skin surfaces if skin integrity is compromised.<sup>1</sup>
- The wearing of masks, eye protection, and face shields in specified circumstances when blood or body fluid exposures are likely to occur is mandated by the OSHA Bloodborne Pathogens Standard.<sup>1</sup>

### COMFORT IS A PRIORITY

- Lightweight design can be worn for long periods of time (< 1.5 oz)
- Adjustable face shield provides flexibility for the worker to move into optimal position
- Includes a nose piece for comfortable and secure fit, even with glasses
- Maximizes airflow to keep the face cooler as compared to blue foam face shields
- Good for sensitive skin and avoiding skin breakouts
- Patented Design

For a no-charge sample of any TIDI product, call **1.800.521.1314** or visit **[www.TIDIproducts.com](http://www.TIDIproducts.com)**

## TIDISHield® Assemble 'n Go™ Face Shield Dispensers Provides protective eyewear when and where you need it.

Designed to maximize airflow and prevent fogging- weighs less than 1.5 ounces! Provides over-the-top, side and front face protection. Disposable, optical grade, anti-static replacement shields. Includes a clear, snap-in nosepiece for support.



### TIDISHIELD® CLEAR SHIELDS FOR ASSEMBLE 'N GO™ FACE SHIELDS (2210-100)

100 Shields  
Use with Visors  
13" x 7"

2210-100



### TIDISHIELD® CLEAR SHIELDS FOR ASSEMBLE 'N GO™ FACE SHIELDS (2210-40)

40 Shields  
Use with Visors  
13" x 7"

2210-40



### TIDISHIELD® ASSEMBLE 'N GO™ FACE SHIELDS (2209)

20 Visors, 20 Nose Pieces and  
40 Shields per Dispenser  
4 Dispensers per Case

2209



### TIDISHIELD® VISORS FOR ASSEMBLE 'N GO™ FACE SHIELDS

(2211-80DSP)  
Assorted Colors  
20 Visors and 20 Nose Pieces  
4 Dispensers per Case

2211-80DSP

### TIDISHIELD® ASSEMBLE 'N GO™ VISORS

(2211-100)  
Assorted Colors  
100 (10 Bags of 10 Visors)



### TIDISHIELD® ASSEMBLE 'N GO™ CLEAR FACE SHIELDS DISPENSER HOLDERS


(7005)  
2 Holders per Case  
8.5" W x 4.75" D X 8.125" H

7005


## ASSEMBLY INSTRUCTIONS


For ultimate under face protection, place visor on the bottom hole, for more room, place visor on the top shield hole.

- 1 Insert nosepiece tab into visor slot, so nosepiece slants away from front of visor. (Nosepiece may be removed to wear over eyeglasses.)

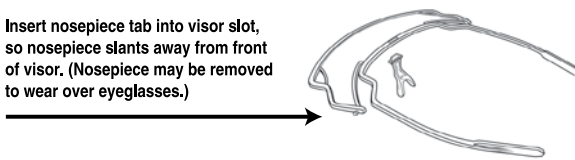
- 2 Starting on one side of the visor, first mount the face shield onto the visor pin (1), then insert the visor's support knuckle (2) through the face shield slot.

Selection of mounting holes provides 3 options for adjusting the face shield angle. Top hole will position the face shield farther from chin, and bottom hole will position it closer to chin.

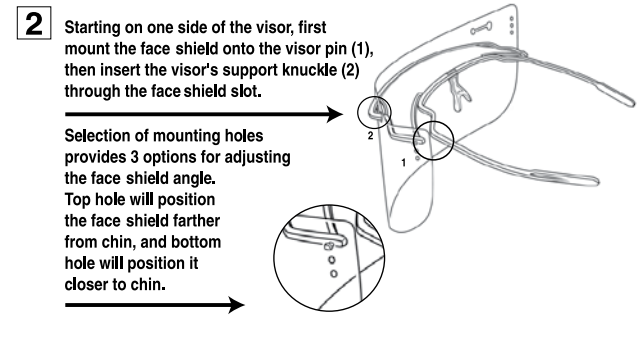
- 3 Repeat step 2 on the opposite side of the visor, pulling the face shield tight and securing it to the opposite support knuckle and pin.

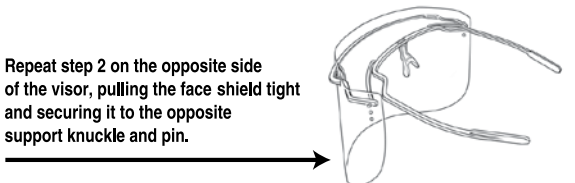



570 ENTERPRISE DRIVE, NEENAH, WI USA 54956  
TOLL FREE 1.800.521.1314  
INTERNATIONAL +1 920.751.4036  
FAX 1.800.837.7770  
WWW.TIDIPRODUCTS.COM



COPYRIGHT © 2016 TIDI PRODUCTS, LLC. ALL RIGHTS RESERVED.