

HC150

The next generation in **CPAP Humidification**



- Ambient Tracking™ to prevent condensation
- DC/AC inverter compatibility
- Improved Design

HC 150 ACCESSORIES		HC 150 SPECIFICATIONS
900HC105	Starter Kit (recommended)	Size 5.4cm x 13.2 cm x 14.6cm (2.1" x 5.2" x 5.7")
HC325D	2x HC325 Chamber	Weight 0.8kg with chamber (1.8 lbs)
HC325T	20x HC325 Chamber	Heater plate range 30°C to 65°C (86°F to 149°F)
900HC007	Mounting Tray	
900HC008	Tube (61cm/24")	
900HC221	Tube (183cm/72")	
HC300	Dishwashable Chamber Kit	

For more information, please contact your local Fisher & Paykel Healthcare representative

Fisher & Paykel have a policy of continuous product improvement and reserve the right to alter specifications without notice. Patented in principal countries of the world.

USA Tel: 1800 446 3908 or +1 949 470 3900 Fax: +1 949 470 3933
INTERNATIONAL P.O.Box 14-348, Panmure, Auckland 6, New Zealand Tel: +64-(0)9-574 0100 Fax: +64-(0)9-574 0158 Email: info@fphcare.com Web Site: www.fphcare.com
AUSTRALIA Tel: +61-(0)3-9879 5022 Fax: +61-(0)3-9879 5232 **FRANCE/BENELUX** Tel: +33-(0)1-64 46 52 01 Fax: +33-(0)1-64 46 52 21
GERMANY/AUSTRIA Tel: +49-(0)7182-93777-0 Fax: +49-(0)7182-93777-99 **UK/IRELAND** (EU Authorised Representative) Tel: +44-(0)1628-626 136 Fax: +44-(0)1628-626 146

Fisher & Paykel
HEALTHCARE

HC150

The next generation in **CPAP Humidification**



Fisher & Paykel
HEALTHCARE

HC150

The next generation in **CPAP Humidification**

Why add Humidification?

Heated humidification improves patient comfort and compliance

Up to 80% of CPAP users complain of upper airway or nasal symptoms. Unfortunately these symptoms may be so severe, patients are unable to continue treatment.

The nasal mucosa is responsible for heating and humidifying the air we breathe. Research has shown the increased airflow of CPAP can overwhelm our natural ability to humidify inspired air. As a result, the nasal mucosa may become dry, irritated and swollen.


The HC150 adds heated humidification to standard CPAP therapy, maintaining mucosal hydration and preventing discomfort.

Ambient Tracking™ to prevent condensation

More humidity, less condensation

During the night, room temperatures often drop. When this occurs, humidification users commonly notice condensation and water droplets forming in the tubing. Many find the associated gurgling noise irritating and are often forced to drain the tubing during the night.

To overcome this, Fisher & Paykel Healthcare has developed the HC150 Humidifier with a unique feature - Ambient Tracking™.

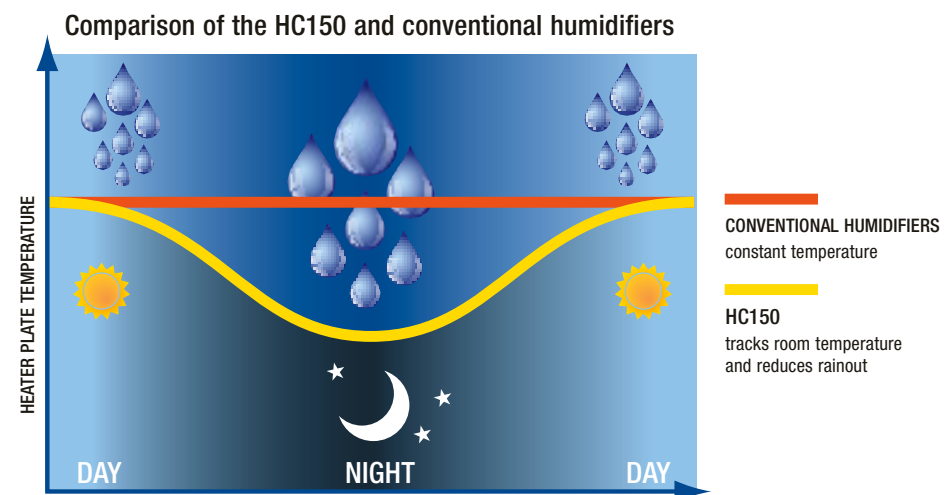

Through advanced electronics, Ambient Tracking™ automatically adjusts the heater plate temperature as the room temperature changes. Condensation is minimized while humidity is maximized.

HC150

CPAP Humidification

Improved Design

- Ambient Tracking™
- Simple humidity controls
- Front mounted on/off switch
- DC/AC inverter compatibility
- Small and compact

Humidification Prevents

- Nasal dryness, congestion and discomfort

Humidification Promotes

- Patient comfort
- Patient compliance



Fisher & Paykel
HEALTHCARE