


# Certofix<sup>®</sup> Quinto

For more safety  
with poly-traumatised patients

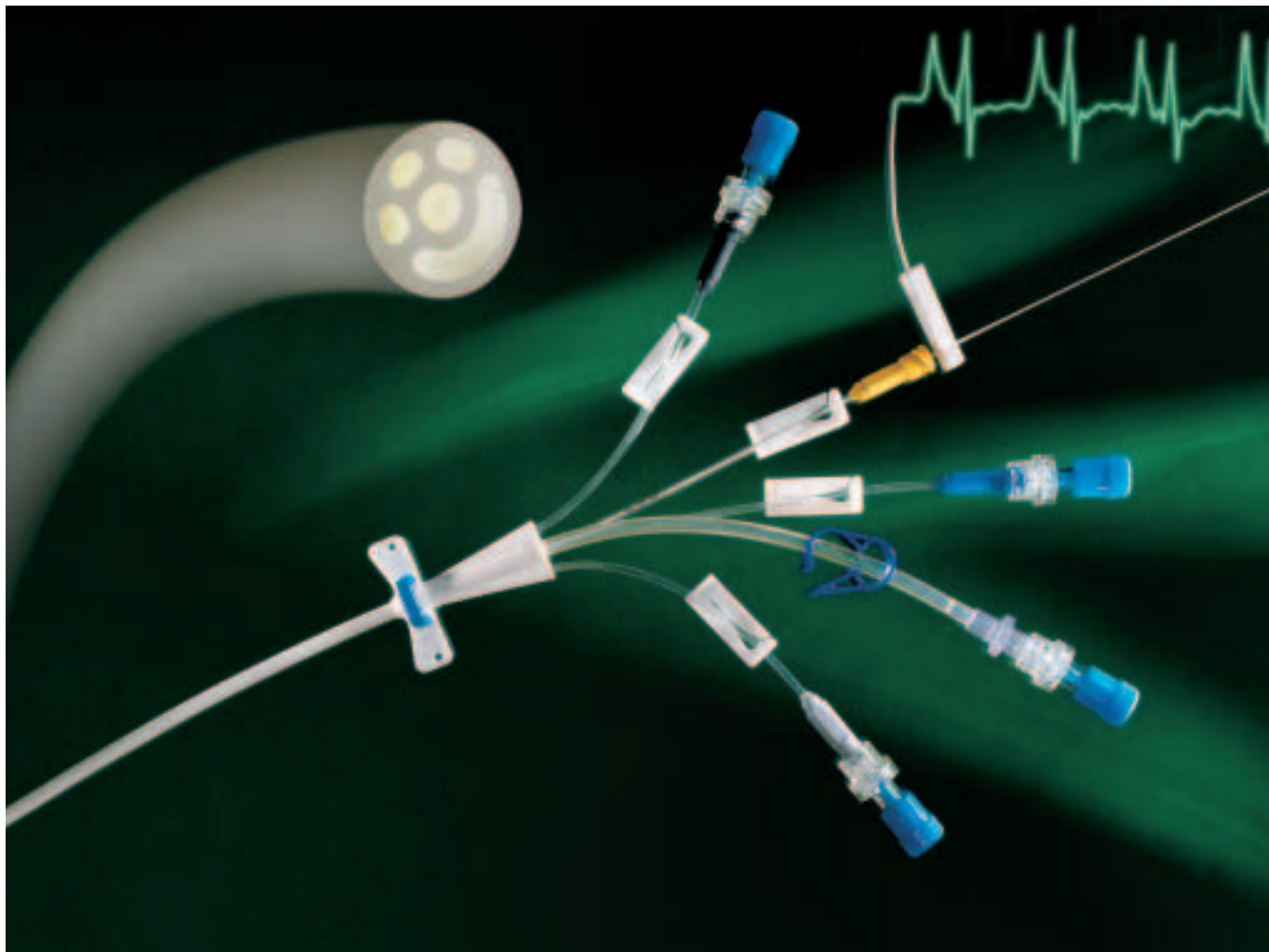

Five-lumen central venous catheters  
for placement using Seldinger technique  
with simultaneous ECG lead option

**B | BRAUN**

# Certofix® Quinto

For more safety  
with poly-traumatised patients

## Five-lumen catheter


- for poly-traumatised patients and for vascular and abdominal surgery
- for an optimal intensive care with complex infusion therapies



- for blood transfusions in case of extensive volume deficiency
- for the avoidance of drug incompatibilities

## Catheter with soft tip

- reduces vein perforations and intima lesions with the flexible, soft polyurethane tip

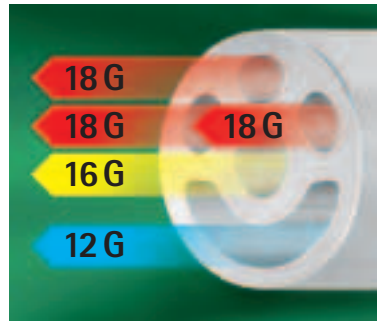


## Optimized section

- 18 G channels for the parallel application of e.g. catecholamines, barbiturates, diuretics etc.
- 16 G channel for the parenteral alimentation
- 12 G channel for transfusions or volume fluid substitution

## Simultaneous positioning check during catheterization via the atrial ECG

- enables immediate corrective measures to be made in case of initial incorrect positioning of the catheter
- helps prevent subsequent complications, particularly the risk of pericardiac perforation and cardiac arrhythmias

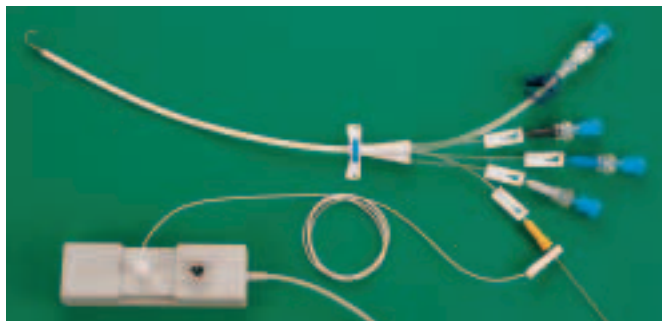

- no X-ray required, no radiation exposure
- no incidents with contrast medium
- effortless identification of the P wave and thus clear verification of position
- no subsequent recatheterization due to incorrect positioning
- reduced costs of equipment and staff

## Product Identification

| Name                     | Code number | Catheter/Lumen Ø G | Length cm |
|--------------------------|-------------|--------------------|-----------|
| Certofix® Quinto SB 1215 | 0416 6833   | 16/18/18/18/12     | 15 cm     |
| Certofix® Quinto V 1215  | 0416 6841   | 16/18/18/18/12     | 15 cm     |
| Certofix® Quinto SB 1220 | 0416 6850   | 16/18/18/18/12     | 20 cm     |
| Certofix® Quinto V 1220  | 0416 6868   | 16/18/18/18/12     | 20 cm     |
| Certofix® Quinto SB 1230 | 0416 6876   | 16/18/18/18/12     | 30 cm     |
| Certofix® Quinto V 1230  | 0416 6884   | 16/18/18/18/12     | 30 cm     |

**B | BRAUN**  
HOSPITAL CARE

B. Braun Melsungen AG  
P.O. Box 11 20  
D-34209 Melsungen

Tel (0 56 61) 71-0  
Fax (0 56 61) 71-45 55  
www.bbraun.com