




Terumo® Dental Needle




Reliable patient comfort

- Ultra-sharp, high-tensile needle
- Reliable fixation
- Bevel indicator
- Individual packaging with clear colour-coding
- Broad selection

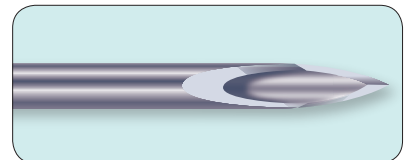

The ultra-sharp lancet point of the Terumo® dental needle, minimises your patient's discomfort and pain. The need for a surface anaesthesia prior to injection can be avoided as well.

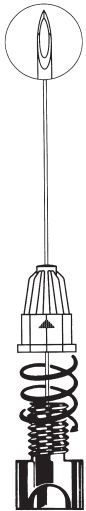

The high-tensile needle is difficult to break, even if the patient should move violently during the injection.

A reliable fixation is guaranteed, as you are able to screw the needle hub directly and securely onto the syringe. No accidental disconnection or leaking will be possible.

A useful bevel indicator, visualised by an arrow on the needle hub, will allow you to align the bevel of your needle for easier insertion.

The individual packaging of each needle in its own plastic sheet, with clear colour coding ensures the sterility of the needle up to the moment of injection. In addition, the colour-coded sealing label provides easy identification of gauge size and facilitates inventory control.


The broad selection of short and long sizes allows you to make the most optimal choice of needle length for conduction anaesthesia as well as infiltration anaesthesia.




Instructions for Use



1. Twist cap to break sealing label and remove the cap.

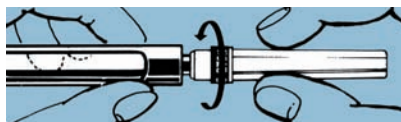

2. Attach needle to syringe by using the cartridge sheath as a wrench to screw the needle hub onto the syringe.

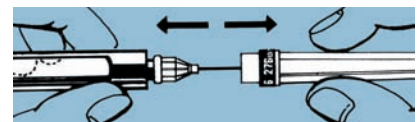

3. Immediately prior to injection, remove the sheath to expose the needle.

Terumo® Dental needle

Gauge	Colour Code	IMPERIAL THREAD (blue hub) Effective needle length							METRIC THREAD (white hub) Effective needle length				
		16 mm	19 mm	21 mm	22 mm	30 mm	35 mm	41 mm	8 mm	13 mm	16 mm	21 mm	35 mm
25G	Red			DN*2521		DN*2530		DN*2541	DN*2508F			DN*2521B	DN*2535B DN*2535F
27G	Yellow			DN*2721	DN*2722	DN*2730	DN*2735	DN*2741			DN*2716F	DN*2721B	DN*2735B
30G	Green	DN*3016	DN*3019	DN*3021	DN*3022					DN*3013B DN*3013F	DN*3016B DN*3016F	DN*3021B DN*3021F	

Unit box: 100 pieces/unit box
Shipping carton: 2000 pieces
Sterilized by ethylene oxide

Needle length at cartridge side: Imperial thread needles: All items = 10.7 mm
Metric thread needles: Items with code B = 13.2 mm
Items with code F = 10.7 mm



Every day, as healthcare professionals, you are working to ensure life and to contribute to the well-being of your patients. Every day, Terumo Europe provides you with the clinical systems, means and solutions to keep life flowing.



TERUMO CORPORATION
44-1, 2-chome
Hatagaya
Shibuya-ku
Tokyo 151-0072
JAPAN
Tel. +81/3 3374.8111
Fax +81/3 3374.8196

TERUMO EUROPE N.V.
Researchpark Zone 2
Haasrode
Interleuvenlaan 40
B-3001 Leuven
BELGIUM
Tel. +32/16.38.12.11
Fax +32/16.40.02.49

TERUMO EUROPE N.V. BENELUX SALES DIVISION
Researchpark Zone 2
Haasrode
Interleuvenlaan 40
B-3001 Leuven
BELGIUM
Tel. +32/16.39.25.80
Fax +32/16.39.25.99
THE NETHERLANDS
Tel. 0800/022.03.96
Fax 0800/022.04.14

LABORATOIRES TERUMO FRANCE S.A.
1, Parc Ariane
Bâtiment Uranus
Rue Hélène Boucher
78284 Guyancourt Cedex
FRANCE
Tel. +33/1 30.96.13.00
Fax +33/1 30.43.60.85

TERUMO DEUTSCHLAND GmbH
Hauptstrasse 87
65760 Eschborn
GERMANY
Tel. +49/61.96.80.230
Fax +49/61.96.80.23.200

TERUMO DEUTSCHLAND GmbH Zweigniederlassung
Wettingen
Bahnhofstrasse 88
CH-5430 Wettingen
SWITZERLAND
Tel. +41/56.419.10.10
Fax +41/56.419.10.11

TERUMO EUROPE ESPAÑA SL
Avda. Juan Carlos I,
N° 13 - 7ª Planta
Edificio Torre La Garena
28806 Alcalá de Henares
(Madrid)
SPAIN
Tel. general:
+34/902.10.12.98
Tel. cust. serv.:
+34/902.20.31.40
Fax general:
+34/902.10.13.58
Fax cust. serv.:
+34/902.10.17.57

TERUMO ITALIA SRL
Via Simone Martini
143/145
I-00142 Roma
ITALY
Tel. +39/06.519.61.420
Fax +39/06.503.04.07

TERUMO SWEDEN AB
Nya Varvet, Byggnad 90
Sven Källfelts Gata 18
SE-426 71 Västra
Frölunda
SWEDEN
Tel. +46/31.748.58.80
Fax +46/31.748.58.90

TERUMO DENMARK Filial af TERUMO SWEDEN AB
Tel. +45/7020.93.80
Fax +45/7020.94.80

TERUMO UK Ltd.
Tamesis
The Causeway
Egham
Surrey
TW20 9AW
UNITED KINGDOM
Tel. +44/1784.476.200
Fax +44/1784.476.201

TERUMO EUROPE N.V. Moscow Representative office
4, Ivana Franko St
121108 Moscow
RUSSIA
Tel. +7/495.234.3303
Fax. +7/495.234.3303

®: Registered Trademark
Published by Terumo Europe N.V.

TERUMO®

We keep life flowing

Addresses are subject to change, please consult our website.

www.terumo.com
www.terumo-europe.com